

COMECTA Eyepiece video camera

ADAPTS DIRECTLY ONTO THE EYEPIECE OF THE MICROSCOPE OR TRINOCULAR HEAD.

Micro video camera 3.1 MP with a USB connection, direct to a PC, includes software.
Part No. **5435030**



COMECTA Light source, cold source, fibre optic "Z-150" CONSTANT COLOR

APPLICATIONS

Microscopy illumination for small samples free from colour alteration.

Especially recommended for microphotography of live biological samples, cultures, micro-organisms, process control, electronic components, dental studies, etc.

FEATURES

Halogen lamp 150 W. Adjustable illumination intensity control. Forced air circulation. Heat absorption filter. Safety transformer fitted. Primary fuse protection.

CONTROL PANEL

Mains switch. Potentiometer controlled light intensity.

MODEL

Part No.	Height / Width / Depth exterior cm	Power W	Weight Kg
5901016	15 19 25	150	4

SPARE Lamp 150 W.
Part No. **5901017**

Supplied with double flexible arm 50 cm long.



COMECTA Microscope accessories and consumables



GLASS SLIDES

Glass slides (smooth): 26 x 76 x 1.5 mm.

Part No. **5600400** Box of 72.

Glass slides (matt end for notes): 26 x 76 x 1.5 mm.

Part No. **5600399** Box of 72.

Cover slip: 18 x 18 mm.


Part No. **5600396** Box of 100.

Cover slip: 22 x 22 mm.

Part No. **5600397** Box of 100.

Cover slip: 24 x 24 mm.

Part No. **5600398** Box of 100.

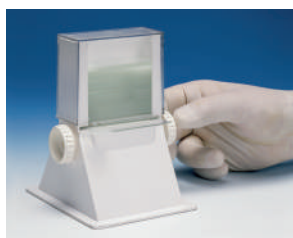

SLIDE DISPENSER

Base made of ABS with a transparent Polycarbonate lid, with rotating manual dispenser of individual slides.

Capacity 50 slides of 76 x 26 mm.

Part No. **1001366**



CLEANING KIT

Suitable for optical microscopes.

Part No. **5313000**

WOODEN MICROSCOPE BOX

Solid wood, microscope box, of 3 mm thickness, treated with varnish, fitted with a metallic door, clasp lock and a carry handle.

Microscope Box suitable for Model M-100FL, XTX-3C and ZTX-20.

420 high x 255 wide x 285 mm deep. Part No. **5313022**
Microscope Box suitable for Model SQF-L, SQF-E, 2000-A, 2000-B and 2000-C.

Internal dimensions: 420 high x 265 wide x 317 mm deep.
Part No. **5313016**

Slide support, Slide box for storing slides, slide staining stands, etc. (See pages 215 to 219).

COMECTA Disk polarimeter "polar-2010" LED

LAURENT TYPE MODEL
LED LAMP ILLUMINATION

APPLICATIONS

Designed for colleges and universities, chemical laboratories, pharmaceutical, clinics, cosmetics and food industry.

FEATURES

Vernier reading definition. Angle $\pm 0,05^\circ$.

Measuring range of optical rotation: $\pm 180^\circ$.

Rotation scale division: 1° .

Measuring lens 4x.

Ocular with dioptres adjustment.

Slide rack compartment length: for tubes up to 220 mm.

Monochromatic lamp wave length: 589,3 nm.

MODEL

Part No	External measurements cm. Height / Width / Depth	Consumption W	Power V / Hz	Weight Kg
5120220	32 43 14	20	110-220 / 50-60	5

SPARE PARTS

Sodium lamp. Part No. **5120221**



NEW
DESIGN