



Blood bank refrigerated cabinets “Blood Bank” A, C and D

FORCED AIR CIRCULATION. FIXED TEMPERATURE +4 °C. STABILITY ±1.5 °C. RESOLUTION 1 °C.
SEVEN DAYS GRAPHIC TEMPERATURE REGISTER.
CAPACITIES FROM 80 TO 480 BLOOD BAGS.



NEW DESIGN

SAFETY:

ALARMS SYSTEMS VISUAL AND ACOUSTIC: LOW TEMPERATURE (+2°C), OVER TEMPERATURE (+6°C), DOOR OPEN, MAINS FAILURE. INDEPENDENT BATTERY BACKUP GRAPHIC TEMPERATURE REGISTER. 7 DAY RECORDER, AUTONOMY 6 MONTHS. ACOUSTIC ALARM ON ALL FAIL CONDITIONS. REMOTE ACOUSTIC ALARM.



Model A



Model C



Model D

APPLICATIONS

AISI 304 storage cabinets specifically designed and manufactured for blood bag storage.

COMMON FEATURES

Model A: interior furniture in painted aluminum and base in stainless steel AISI 304, exterior in white enameled steel.

Models C and D: Stainless steel. AISI 304. external cabinet and interior lining.

Double glazed antimist, reversible lockable door, easy to change, self close door if left ajar. Hermetically sealed compressor mounted on shock absorbers with evaporator and forced air condenser.


Fan air circulation maintaining homogeneous temperature distribution.

Temperature sensors are located and calibrated to read simulated blood bag conditions.

Independently powered circular graphic register.

Microprocessor controlled temperature control via a display that can configure the following parameters:

- Adjust SET.
- Adjust alarm.
- Calibration.

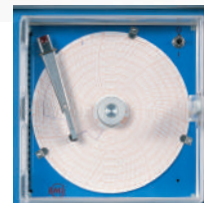

Graphic recorder

OPTIONAL COMPLEMENTS

ETHERNET thermometer module

Temperature Record, with control of programmable warnings and graphic representation. It must be installed at the factory.

ETHERNET connection is needed at the place of use.

Part No. **2101627**

Data logger with digital display

PT100 type probe with two meters of cable. Measuring range of -100 °C. at +200 °C. With GSM module for data transmission.

Part No. **2101520**

SPARE PARTS

For 2101519. Part No. **1001864** Additional tray.

For 2101512. Part No. **1001862** Additional tray.

For 2101514. Part No. **1001863** Set 2 trays + additional central guide

MODELS	Part No.	Capacity litres	Capacity bags	Height / Width / Depth (interior) cm			Height / Width / Depth (exterior) cm			Power W	Power HP	Number of shelves	Weight Kg
Blood Bank A	2101519	235	80	109	51	55	156	61	62	300	1/7	4	90
Blood Bank C	2101512	560	200	153	54	66	213	80	79	497	3/8	4	165
Blood Bank D	2101514	1240	480	153	124	68	213	139	79	1034	5/8	8	212

They are supplied with trays and wheels.