


Dry bath

NEW

compact design, 1-block-model, up to 150°C, various blocks available

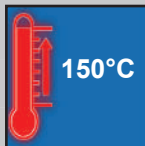

H-20 up to 150°C

Ideal for:

- enzyme reactions, inactivation of sera, incubation and other laboratory procedures
- multifunctional dry block heater

Features:

- new design: lighter and more compact than the previous model
- user's self-compensation function to offset between real temperature and indicated temperature ($\pm 10.0^\circ\text{C}$)
- ergonomic, compact design & modular block system
- internal PPS bath ensures best insulation, block made from anodized aluminium
- precise temperature control
- digital timer with delay time and operation time
- convenient integrated timer for time-sensitive incubations
- block removal tool with screw thread included
- storage function for the set values of temperature and timer
- various block types (optional)
- backlight LCD display
- CE certified and unique serial number for tracing

Safety mechanism:

- locking mode
- overheat and over-current protection
- sensor error detection

Controller:

- digital fuzzy controller with **Jog-Dial** and **main touch button**




H-20 set with block and block removal tool

Model	H-20	
Block	interchangeable block (90 x 120 x 60 mm) for various tubes	
Standard block	-	2x 48 holes for 1.5/2 ml centrifuge tubes
Temperature range and Accuracy	ambient to $+5^\circ\text{C}$ to 150°C , ± 0.1 at 37°C , $\pm 0.2^\circ\text{C}$ at 90°C , $\pm 0.4^\circ\text{C}$ at 135°C	
Heat up speed & sensor	approx. $5^\circ\text{C}/\text{min}$, PT100	
Timer with alarm	99hr 59min, alarm (error status and timer end)	
Dimension	unit: 178 x 226 x 138 mm, block: 90 x 120 x 60 mm	
Packing size	230 x 290 x 245 mm, 2.2 kg	230 x 355 x 245 mm, 3.5 kg
Power supply	1 Phase AC 120V, 60 Hz or AC 230V, 50/60Hz	
Power consumption	300 W	
Order number 230V	Unit: DH.WHB02000	Set: DH.WHB02100
Order number 120V	Unit: DH.WHB02001	Set: DH.WHB02101

Set includes: Unit
Block
Block removal tool

H-20
BLC548 for 48x 1.5/2 ml tubes

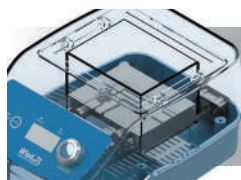

High quality molded heater ensures reliable and consistent temperature




Maximum heat insulation and highest chemical resistance by PPS materials

Accessories: Dry bath blocks and removal tool on page 33

Dry bath

two blocks, individual controller for each block, up to 150°C, various blocks available

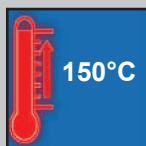

HB-96D 2 blocks up to 150°C

Ideal for:

- enzyme reactions, inactivation of sera, incubation and other laboratory procedures
- multifunctional dry block heater

Features:

- user's self-compensation function to offset between real temperature and indicated temperature ($\pm 10.0^{\circ}\text{C}$)
- unit made of stainless steel and powder-coated steel, block made of anodized aluminium
- digital timer function
- block extraction tongs "BLE101" & block cover are included
- storage function for the set values of temperature and timer
- alarm function: error status and timer end
- various kinds of blocks are offered (optional)
- dual block and 2 molded heater (300W)
- backlight LCD display
- CE certified and unique serial number for tracing

Safety mechanism:

- locking mode
- overheat and over-current protection
- sensor error detection

Controller:

- dual digital fuzzy controller with *Jog-Shuttle-switch*



HB-96D Set with blocks and tong

Model	HB-96D	
Block	interchangeable block (90 x 120 x 60 mm) for various tubes	
Standard block	-	2x 48 holes for 1.5/2 ml centrifuge tubes
Temperature range and Accuracy	ambient to $+5^{\circ}\text{C}$ to 150°C , ± 0.1 at 37°C , $\pm 0.2^{\circ}\text{C}$ at 90°C , $\pm 0.4^{\circ}\text{C}$ at 135°C	
Heat up speed & sensor	approx. $5^{\circ}\text{C}/\text{min}$ PT 100	
Timer with alarm	99hr 59min, alarm (error status and timer end)	
Dimension	unit: 285 x 250 x 155 mm, 2 x block: 90 x 120 x 60 mm	
Packing size	400 x 270 x 200 mm, 6 kg	
Power supply	1 Phase AC 120V, 60 Hz or AC 230V, 50/60Hz	
Power consumption	600 W	
Order number 230V	Unit: DH.WHB00198	Set: DH.WHB00199
Order number 120V	Unit: DH.WHB02198	Set: DH.WHB02199

Set includes: Unit HB-96D
Block 2x BLC548 for 48x 1.5/2 ml tubes
Tong BLE101

Accessories: Dry bath blocks and tong on page 33

Blocks for dry baths

- anodized aluminium
- 90 x 120 x h60 mm



For	0.2 ml PCR tubes	0.5 ml PCR tubes	1.5/2 ml tubes	Ø12 mm tubes	Ø13 mm tubes
Holes	96 holes Ø: 6 mm depth: 15.5 mm conical bottom	63 holes Ø: 7.7 mm depth: 29.5 mm conical bottom	48 holes Ø: 11 mm depth: 36 mm conical bottom	48 holes Ø: 12 mm depth: 35 mm round bottom	30 holes Ø: 13 mm depth: 47 mm round bottom
Article	BLC596	BLC570	BLC548	BLT248	BLT335
Order number	DH.WHB510596	DH.WHB510570	DH.WHB511548	DH.WHB521248	DH.WHB521335



For	Ø15 mm tubes	Ø16.5 mm tubes	Ø18 mm tubes	15 ml conical tubes	50 ml conical tubes
Holes	30 holes Ø: 15 mm depth: 47 mm round bottom	24 holes Ø: 16.5 mm depth: 47 mm round bottom	24 holes Ø: 18 mm depth: 47 mm round bottom	15 holes Ø: 16 mm depth: 47 mm conical bottom	8 holes Ø: 28.7 mm depth: 47 mm conical bottom
Article	BLT535	BLT628	BLT828	BLC515	BLC008
Order number	DH.WHB521535	DH.WHB521628	DH.WHB521828	DH.WHB531515	DH.WHB535008

Block removal tool with screw thread, for adding and removing blocks

Order number

DH.WHB500102

NEW



Tong, for adding and removing blocks

Article

BLE101

Order number

DH.WHB500101

