


Bath | water bath

extraction water bath, 4 / 6 / 8 places, up to 100°C ±0.2°C, holding frame and covers included

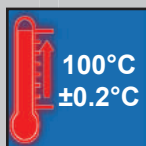

100°C
±0.2°C

WEB-6 with holding frame (included) and concentric ring covers (included)



WEB-4 with concentric ring covers (included)

WEB 100°C, 4/6/8 places, digital

Ideal for:

- gentle extraction, soxhlet extraction and determination of COD

Features:

- available with 4,6 or 8 extraction places
- user's self-compensation function to control difference between real temperature and indicated temperature (±10.0°C)
- concentric ring cover set Ø 110 mm and stainless steel holding frame included, drain valve mounted
- stainless steel interior and powder-coated steel body
- storage function for the set values of temperature and timer
- backlight LCD display
- CE certified and unique serial number for tracing

Safety mechanism:

- overheat and over-current protection
- sensor error detection
- locking mode
- low water level warning and automatic stop function

Controller:

- digital fuzzy controller with **Jog-Shuttle-switch** (turn + push)

Model	WEB-4	WEB-6	WEB-8
Temperature range & accuracy	room temperature +5°C - 100°C, ±0.2°C		
Temperature uniformity & sensor	±1.0°C, PT100		
No. of extraction holes	4 places	6 places	8 places
Heating power	1.4 kW	1.8 kW	2 x 1.4 kW
Timer & alarm	99hr 59 min (delay & continuous run), error status & timer end		
Internal dimensions (W x D x H)	610 x 160 x 160 mm	910 x 140 x 160 mm	610 x 290 x 160 mm
External dimensions (W x D x H)	720 x 245 x 275 mm	975 x 245 x 275 mm	720 x 375 x 275 mm
Packing size & gross weight	710 x 310 x 300 mm, 18 kg	1,030 x 320 x 370 mm, 25 kg	880 x 420 x 300 mm, 34 kg
Power supply	1 Phase AC 120V, 60 Hz or AC 230V, 50/60 Hz		
Order number 230V	DH.WEB01004	DH.WEB01006	DH.WEB01008
Order number 120V	DH.WEB02004	DH.WEB02006	DH.WEB02008

Accessories:



Holding frame set, stainless steel

Measurements	Article	Order number
for WEB-4, 2x800mm pole, 2x700mm pole	EBF100	DH.WEB21004
for WEB-6, 2x800mm pole, 2x1.020mm pole	EBF200	DH.WEB21006
for WEB-8, 4x800mm pole, 4x740mm pole	EBF300	DH.WEB21008

Scope of supply:

- WEB water bath
- holding frame
- concentric ring covers



Lid with concentric ring, stainless steel

	Article	Order number
for all WEB models	EBC100	DH.WEB11000