


Bath | oil bath

oil bath, digital, 6 / 11 / 22 l, up to 250°C ±1.5°C, flat lid included

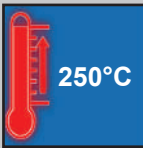

WHB oil bath, 250°C, digital

Ideal for:

- oil bath for general purpose with high temperature

Features:

- user's self-compensation function to control difference between real temperature and indicated temperature
- stainless steel bath & lid for superior durability & high thermal efficiency / consistency
- digital timer function (99hr 59min): delay time and operation time can be controlled
- storage function for the set values of temperature and timer
- backlight LCD display
- CE-certified and unique serial number for tracing

Safety mechanism:

- locking mode
- over-temperature and over-current protection
- sensor error protection

Controller:

- digital fuzzy controller with **Jog-Shuttle-switch** (turn + push)



WHB-11

Flat lid (included)

Model	WHB-6	WHB-11	WHB-22
Capacity	6 l	11 l	22 l
Temperature range & accuracy	room temperature +5°C - 250°C, ±1.5°C		
Temperature resolution	0.1°C-display, 0.1°C-control		
Temperature uniformity & probe	±3.0°C, PT100		
Heating power	1.5 kW	2 kW	2.5 kW
Timer & alarm	99hr 59 min (delay & continuous run), error status & timer end		
Drain valve	mounted		
Effective space	270 x 170 x 110 mm	300 x 250 x 110 mm	400 x 300 x 140 mm
Internal dimensions (W x D x H)	270 x 180 x 150 mm	300 x 250 x 150 mm	400 x 300 x 180 mm
External dimensions (W x D x H)	425 x 305 x 340 mm	425 x 375 x 340 mm	525 x 425 x 370 mm
Packing size & gross weight	470 x 400 x 500 mm, 23 kg	600 x 490 x 500 mm, 25 kg	600 x 490 x 500 mm, 39 kg
Power supply	1 Phase AC 120V, 60 Hz or AC 230V, 50/60 Hz		
Order number 230V	DH.WHB00106	DH.WHB00111	DH.WHB00122
Order number 120V	DH.WHB00206	DH.WHB00211	DH.WHB00222